




PHOENIX

Finds every leak. Now with an even finer nose.

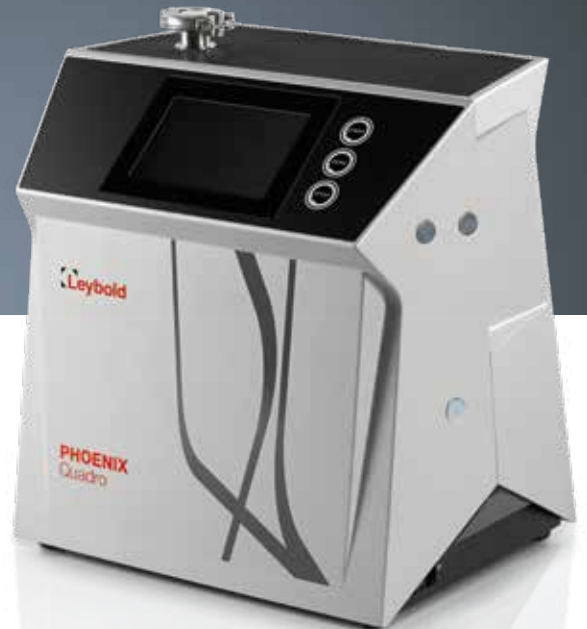

The new generation of leak detectors. Where precision meets speed, with your convenience in mind.

Often the success of any business depends on quality, and having the right tools and equipment is critical to ensuring top performance and output across R&D departments, as well as in fast-paced production and manufacturing facilities. We also see applications becoming more specialized, processes being more demanding, yet time for analysis and testing becoming less.

However, when it comes to vacuum system leak detection you can rely on Leybold and the newest generation of PHOENIX leak detectors, where precision meets speed, with your convenience in mind.

A new degree of precision

Whether you are in a research lab, or production of high-class goods, quality depends on how accurately you detect deviations. Now, help is at hand: the new PHOENIX generation delivers the highest helium sensitivity for unbeatable reliability and stability.

A new level of speed

If you are looking to accelerate the leak test of your system or speed up quality control in manufacturing, time is of the essence. The new PHOENIX leak detectors feature an improved measurement system with faster response times than ever before.

A new dimension of comfort and choice

The new PHOENIX offers a choice of three models, so you are sure to find the right leak detector for your needs, with simple operation, integrated touch-screen and Wi-Fi options. New software functions let you add an external display to monitor and run measurements from anywhere for added flexibility. And thanks to internal and external data logging, no measured value is ever lost.



“The new PHOENIX leak detector really wins you over with its outstanding accuracy, great performance and ease of use. It’s a leak detector that you enjoy working with – because it’s intuitive and gets you every leak.”

Dr Sina Forster, Product Manager Leybold



A touch easier

Say goodbye to complicated menus on illegible displays. The new PHOENIX has state-of-the-art multicolored touchscreens with an intuitive menu structure. Hello, ease of use!



No fuss: easy login

The new PHOENIX models are the only portable detectors on the market featuring easy login via Wi-Fi. Simply use your smartphone or tablet – no software or app to download.

Extremely sensitive

All PHOENIX leak detectors feature an improved measurement system, with a new ULTRA mode for excellent helium sensitivity. Even the smallest leaks can be detected rapidly. The result: faster response times than ever before.





Dozens of applications. Three strong solutions.

Whoever works with vacuum needs a leak detector. Ideally one that delivers precise results, time after time. One that offers easy data logging as well as remote operation. One that stands out by fitting in, seamlessly, to your system.

Whatever your application, our new PHOENIX series has you covered. Discover all the advantages of the latest generation!

Research & Development

- Universities and institutes
- Research centers
- Labs
- R&D equipment manufacturing

- Compact systems enable optimum use of available space
- High mobility and flexibility
- Easy and intuitive operation
- Intelligent reminders and convenient data logging options
- Change easily between different users: Set up user profiles and log in quickly
- Change easily between different experiments: Create the test set-up and save up to 10 parameter sets
- Easy data analysis and documentation using a variety of interfaces

Your best bet: PHOENIX Quadro or Vario on a mobile cart

UHV Applications

- Particle accelerators
- UHV components and systems
- Space simulation

- Smallest detectable leak rate of 5×10^{-12} mbar·l/s for all models
- Highest accuracy and reliability of measurements
- Very low background and intelligent zero function for quick measurements of low leak rates
- High pumping speed for fast measurements when it comes to large volumes
- Dry versions meet the need for oil-free systems and processes

PHOENIX Vario or Magno dry are excellent choices



The PHOENIX offers all the necessary connections for accessories including the sniffer line or the partial flow system. Furthermore, it has easily accessible USB ports and additional controlling interfaces.



We offer a range of interface modules that offer suitable communication standards for different needs. Choose between an I/O module, Profibus®, Profinet®, Devicenet® or Ethernet/IP.



The PHOENIX line has an improved service and maintenance concept. The filter can be changed from the outside, without removing the cover. This lets you save time and focus on your measurements!



Industrial Applications

- Industrial systems and machinery
- Metallurgy
- Industrial furnaces
- Electrical industry
- Highest pumping speed is perfectly suited for large volumes
- Harsh processes make robustness relevant: Proven industrial design with rugged components stands up to heavy use
- Big systems require remote control options: Use the RC 310 or the Wi-Fi function and log in with a smartphone or tablet
- Fast operation: Quick run-up and ready to start within seconds

Your perfect fit: PHOENIX Vario

Industrial production

- Analytical instruments
- Medical devices
- Food and packaging
- Automotive and aviation
- Meet high quality expectations: Reliable measurements for consistent and highly accurate test results
- Ensure the highest level of safety: Reliable pass-or-fail decisions
- Greater efficiency: Save parameter sets for different test parts and upload them quickly and easy
- Convenient communication options and various interfaces
- Highest helium pumping speed for fast cycling

PHOENIX Quadro dry is the ideal “nose” here

Semiconductor

- Semiconductor fabrication
- Flat panel display
- Time-saving measurements: Quick to start up, ready to go within seconds
- Small systems offer highest mobility in large facilities
- Outstanding accuracy and repeatability for reliable leak rate measurements

Choose PHOENIX Quadro dry or PHOENIX Vario

Find your perfect PHOENIX



To become a top dog, you have to beat the rest of the pack. Meet our “noses” and their unbeatable features. Then decide which one suits your application best.

The new PHOENIX Quadro and Quadro dry

- Available with oil-sealed or dry backing pump
- Detection sensitivity improved by turbo pump's optimized helium pumping speed
- Pumping speed during evacuation from 2.5 m³/h up to 3 m³/h
- More data flexibility: USB and bus interface options

Technical data	PHOENIX Quadro	PHOENIX Quadro dry
Smallest detectable helium leak rate		
Vacuum operation	mbar · l · s ⁻¹	≤ 5 · 10 ⁻¹²
Sniffer mode	mbar · l · s ⁻¹	< 1 · 10 ^{-9 1)}
Smallest detectable hydrogen leak rate		
Vacuum operation	mbar · l · s ⁻¹	≤ 1 · 10 ⁻⁸
Sniffer mode	mbar · l · s ⁻¹	< 1 · 10 ⁻⁷
Units of measurement (selectable)		
Pressure		mbar, Pa, atm, Torr
Leak rate		mbar · l · s ⁻¹ , Pa · m ³ · s ⁻¹ , Torr · l · s ⁻¹ , atm · cc · sec ⁻¹ , sft ³ /yr
Sniffer mode		ppm, g/a eq, oz/yr eq
Leak rate display area	mbar · l · s ⁻¹	1 · 10 ⁻¹² to 1 · 10 ⁻¹
Max. inlet pressure	mbar	15
with partial flow pump set	mbar	1000
Pumping speed during evacuation process		
50 Hz	m ³ · h ⁻¹	2.5
60 Hz	m ³ · h ⁻¹	3.0
Helium pumping speed in vacuum mode		
GROSS Modus	l/s	0.4
FINE Modus	l/s	1.2
ULTRA Modus	l/s	> 3.1
Time constant of the leak signal (blank flanged, 63% of the final value)	s	< 1
Run-up time (after switching on)	s	≤ 110
Mass spectrometer		180° magnetic sector field
Ion source		2 yttrium/iridium long-term cathodes
Detectable masses	amu	2, 3 and 4
Test port		
Inlet flange	DN	25 ISO-KF
Dimensions	(W x H x D) mm	495 x 475 x 318
Weight	kg	41
Available languages		English, German, Chinese, Japanese (Katakana), Korean, Russian, French, Italian, Spanish and Polish

Order data	PHOENIX Quadro	PHOENIX Quadro dry
	Part no.	Part no.
PHOENIX Quadro EU 230 V, 50/60 Hz	250000V02	-
PHOENIX Quadro US 115 V, 60 Hz	251000V02	-
PHOENIX Quadro JP 100 V, 60 Hz	251100V02	-
PHOENIX Quadro 100 – 240 V, 50/60 Hz	-	250001V02
PHOENIX Magno EU 200 – 230 V, 50/60 Hz		
PHOENIX Magno JP/US 100 – 115 V, 50/60 Hz		
PHOENIX Vario		

¹⁾ Depending on the used sniffer line and ambient conditions. Leak rate display after zero: 1·10⁻⁸ mbar·l/s when using the SL 300 or SL 301 Sniffer line.



Meet the all-new PHOENIX series.

Your benefits at a glance:

- Highest accuracy and reliability thanks to proven measurement system
- Fastest helium pumping speed allows most efficient measurements
- Consistency and convenience: Save all parameters and copy them to other PHOENIX leak detectors
- Innovative touchscreen with intuitive menu structure
- Handy and portable: very small dimensions and low weight
- Log in using Wi-Fi on your smartphone or tablet. No app or software needed!
- Easy internal and external data logging: no measured values are ever lost
- HDMI interface lets you check leaks from anywhere by adding external displays
- Find the perfect fit: three different models to better suit your needs



Leybold GmbH
Bonner Str. 498 · D-50968 Köln

T +49 (0) 221-347-0
F +49 (0) 221-347-1250
info@leybold.com

www.leybold.com
www.phoenix-series.com

